



**Santec surgical sutures
for cardiovascular and general surgery**

**Santec medical products
Production and distribution of single use medical
devices for infusion, transfusion, intensive care
and anesthesia applications.**

CE 1275

**INDIVIDUAL
ORDERS ON
DEMAND**

**san
tec**

**san
tec**

Production - Distribution - Product Management

Santec Medicalprodukte GmbH
Mainstrasse 32
63868 Großwallstadt
Germany
Telephone +49(0)6022 508951-0
Telefax +49(0)6022 508951-29
info@santec-medical.de
www.santec-medical.de

Detlev-Sven Neumann, General Manager
Wolfgang Leichtenstern, General Manager,
Sales Southern Germany

Order Handling

Bernadette Wießmann
Telephone +49(0)6022 508951-0
Telefax +49(0)6092 508951-29
info@santec-medical.de

Sales Eastern Germany

Brigitte Behr-Ostheimer
Mobil +49(0)160 9099 8302
brigitte.behr-ostheimer@santec-medical.de

Sales North Rhine-Westphalia


Horst Derksen
Mobil +49(0)177 1652 029
horst.derksen@santec-medical.de

Sales Northern Germany

Uwe Giebeler
Mobil +49(0)151 5889 0465
uwe.giebeler@santec-medical.de

For foreign clients please refer to our headquarters.

implant Chordae Loop

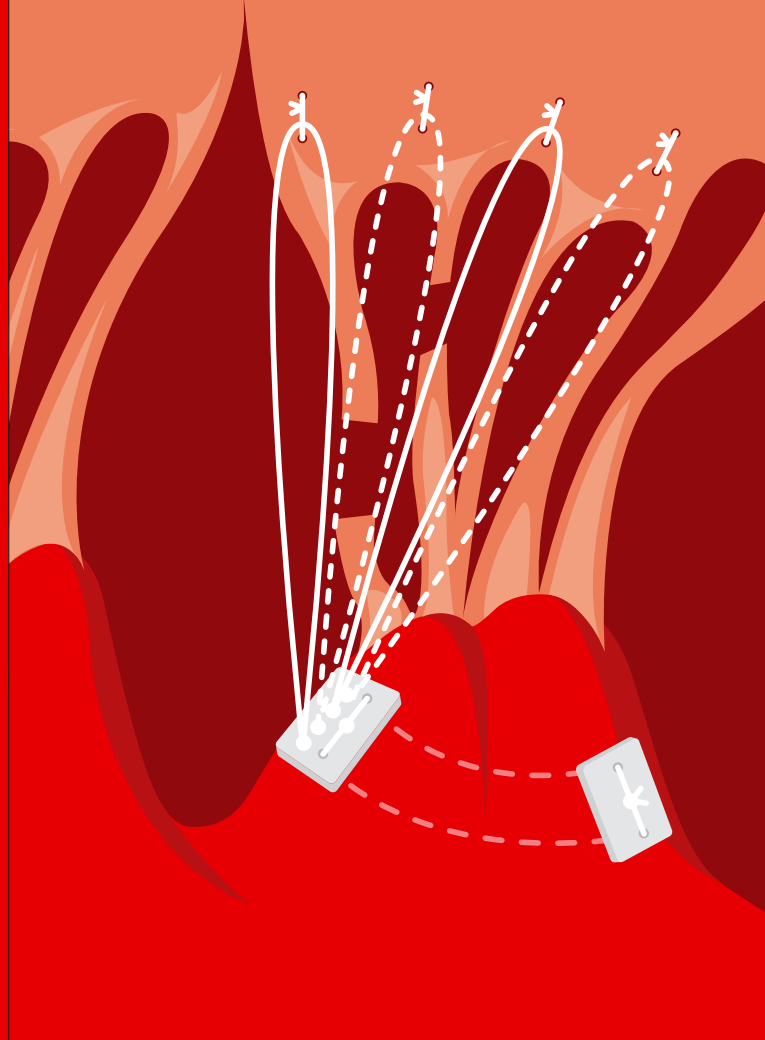

santec medicalprodukte gmbh

www.santec-medical.de

santec medicalprodukte gmbh



implant *Chordae Loop*

The Chordae Loop is a non absorbable implant for Chordae repair/replacement used for mitral valve reconstruction. The implant is based on an end-to-end ePTFE suture, size CV-4 with swaged needles PH17 on both ends. The implant consists of 4 identical loops knotted on the upper surface of a PTFE pledget, the suture ends protrude on the lower surface of the pledget.

The pledget with the knotted loops is sewn onto the papillary muscle and a second pledget serves as a suture support structure. The loops are to be attached to the free margin of the leaflet.

Product Specifications:


- ePTFE suture CV-4 implant with 4 loops of identical length
- Both suture ends are swaged with a PH17 needle
- The sutures are attached to a pledget "P" measuring (6x3 mm)
- A second pledget "P" is provided as a support structure for the suture
- Total suture lengths vary from 35 - 42 cm, according to the respective loop lengths
- Available loop size length (in mm):

10 = REF 710	20 = REF 703
12 = REF 711	22 = REF 704
14 = REF 700	24 = REF 712
16 = REF 701	26 = REF 713
18 = REF 702	

Packaging options:

- 4x1 implant, wrapped individually, sterile, semi box
- 12 x 1 implant on demand

Pledgets - Suture Support Structures

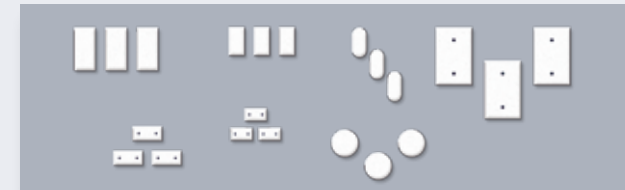

Santec supplies non absorbable pledgets in various sizes and qualities for cardiovascular and general surgery.

The pledgets are made from PTFE-felt (Teflon) and are certified for use as medical implants and as support for medical sutures in delicate tissues.

Surgical Sutures

Santec supplies a variety of surgical sutures, including absorbable and non absorbable sutures in combination with various surgical needle types. Santec can configure and modify the products precisely according to our clients' specifications and needs.

The customers define the combination of suture, needle and pledget.